

Rapid fFN® Cassette Kit

REF 01200

For *In Vitro* Diagnostic Use Only
Store at room temperature (15° to 30°C / 59° to 86°F).

Rx only

For Professional Use Only

INTENDED USE

The Rapid fFN Cassette for use in the TLiQ® System is an *in vitro* diagnostic device for the detection of fetal fibronectin in cervicovaginal secretions to be used as an aid to rapidly assess the risk of preterm delivery in ≤ 7 or ≤ 14 days from the time of cervicovaginal sample collection in pregnant women with signs and symptoms of early preterm labor, intact amniotic membranes and minimal cervical dilatation (< 3 cm), sampled between 24 weeks, 0 days and 34 weeks, 6 days of gestation.

The Rapid fFN test is further indicated for use in conjunction with other clinical information as an aid to rapidly assess the risk of preterm delivery in ≤ 34 weeks, 6 days when a cervicovaginal sample is obtained during a routine prenatal visit between 22 weeks, 0 days and 30 weeks, 6 days of gestation in women with a singleton pregnancy.

The Rapid fFN test represents a significant and critically needed improvement in the ability to manage preterm labor that may result in preterm delivery.

CONTRAINDICATIONS

The Rapid fFN test **should not** be used for symptomatic women with one or more of the following conditions:

- advanced cervical dilatation (≥ 3 centimeters)
- rupture of amniotic membranes
- cervical cerclage
- moderate or gross vaginal bleeding

Delivery typically occurs imminently when the cervix is dilated more than 3 centimeters or if the amniotic membranes are ruptured. Additional diagnostic testing is usually not necessary to confirm risk for women with advanced cervical dilatation or rupture of amniotic membranes. Moderate or gross vaginal bleeding is an independent risk factor for preterm delivery and may be associated with other severe obstetrical or medical problems. Clinical attention should be focused on identification of the origin of bleeding rather than immediate assessment of delivery risk. At this time, information is insufficient regarding the association of vaginal fetal fibronectin expression to delivery for women with cervical cerclage.

The Rapid fFN test **should not** be used for asymptomatic women with one or more of the following conditions:

- multiple gestations, e.g., twins
- cervical cerclage
- placenta previa (partial or complete)
- sexual intercourse in the preceding 24 hours

At this time, information is insufficient regarding the association of cervicovaginal fetal fibronectin expression to delivery for asymptomatic women with HIV/AIDS, multiple gestations or cervical cerclage.

RapidfFN
for the TLiQ® System

Rapid fFN® Cassette Kit

SUMMARY AND EXPLANATION OF THE TEST

Of the approximately 4,000,000 deliveries that occur annually in the United States, about 400,000 are premature. Preterm delivery, defined by the American College of Obstetricians and Gynecologists as delivery prior to the 37th week of gestation, is responsible for the majority of non-chromosomal perinatal morbidity and mortality (1–4). Symptoms of threatened preterm delivery include uterine contractions, change of vaginal discharge, vaginal bleeding, backache, abdominal discomfort, pelvic pressure, and cramping. Diagnostic modalities for identification of threatened preterm delivery include uterine activity monitoring and performance of a digital cervical examination, which allows estimation of cervical dimensions. These methods have been shown to be limited, as minimal cervical dilatation (< 3 centimeters) and uterine activity occur normally and are not necessarily diagnostic of imminent preterm delivery (5,11,13). While several serum biochemical markers have been evaluated, none have been widely accepted for practical clinical use (6,7,20).

Fetal fibronectin (fFN), an isoform of fibronectin, is a complex adhesive glycoprotein with a molecular weight of approximately 500,000 daltons (8,9). Matsuura and co-workers have described a monoclonal antibody called FDC-6, which specifically recognizes III-CS, the region defining the fetal isoform of fibronectin (8,9). Immunohistochemical studies of placentae have shown that fFN is confined to the extracellular matrix of the region defining the junction of the maternal and fetal units within the uterus (5,10).

Fetal fibronectin can be detected in cervicovaginal secretions of women throughout pregnancy by use of a monoclonal antibody-based immunoassay. Fetal fibronectin is elevated in cervicovaginal secretions during the first 24 weeks of pregnancy but is diminished between 24 and 34 weeks in normal pregnancies. The significance of its presence in the vagina during the first 24 weeks of pregnancy is not understood. However, it may simply reflect the normal growth of the extravillous trophoblast population and the placenta. Detection of fFN in cervicovaginal secretions between 24 and 34 completed weeks gestation is reported to be associated with preterm delivery in symptomatic (5,11-15) and asymptomatic pregnant women (16–19).

PRINCIPLE OF THE TEST

The Rapid fFN Cassette is a lateral flow, solid-phase immuno-chromatographic assay. The cervicovaginal specimen is extracted into a buffer and a 200 µL sample is dispensed into the sample application well of the Rapid fFN Cassette. The sample flows from an absorbent pad across a nitrocellulose membrane via capillary action through a reaction zone containing murine monoclonal anti-fetal fibronectin antibody conjugated to blue microspheres (conjugate). The conjugate, embedded in the membrane, is mobilized by the flow of the sample. The sample then flows through a zone containing goat polyclonal anti-human fibronectin antibody which captures the fibronectin-conjugate complexes. The remaining sample flows through a zone containing goat polyclonal anti-mouse IgG antibody which captures unbound conjugate, resulting in a control line. After 20 minutes of reaction time, the intensities of the test line and control line are interpreted with the TLiQ® Analyzer.

PRECAUTIONS AND WARNINGS

Note: *Transport specimens at 2° to 25°C, or frozen. Specimens are stable for up to eight (8) hours at room temperature. Specimens not tested within eight hours of collection must be stored refrigerated at 2° to 8°C and assayed within three (3) days of collection, or frozen and assayed within three (3) months to avoid degradation of the analyte. Specimens arriving frozen may be tested as described below (subject to a single freeze-thaw cycle only).*

1. For *in vitro* diagnostic use only.
2. Carefully follow the instructions and procedures described in this insert.
3. **Test results may not be interpreted visually and must be based on the use of the TLiQ System.**
4. **Do not use glass tubes or glass pipettes, as fetal fibronectin binds to glass. Tubes and pipettes of polypropylene or polyethylene are acceptable.**
5. Do not mix materials from different kit lots.
6. Do not use cassettes or controls past their expiration dates.
7. Do not use controls if they are cloudy or discolored. Avoid cross-contamination of reagents. Use a new pipette tip for each control or patient sample. Recap controls tightly with the correct color-coded caps.
8. Handle cassettes with care; do not touch, scratch, or compress membrane materials in the Rapid fFN Cassette.
9. Source material used to prepare the controls is of human origin. The donors were tested and found to be negative for HIV 1, HIV 2, and HCV antibody and hepatitis B surface antigen (HBsAg) using established methods. No known test method can offer total assurance that HIV, hepatitis C virus, hepatitis B virus, or other infectious agents are absent. **Handle the controls and all patient specimens as if potentially infectious.**

*Rapid*fFN®
for the TLiQ® System

Rapid fFN® Cassette Kit

10. Labels (e.g., bar code labels) can be placed on the thumb grip area of the cassette. Do not place labels on an area of the cassette that will be inserted into the TLiQ® Analyzer.

STORAGE

The Rapid fFN Cassettes should be stored at room temperature (15° to 30°C / 59° to 86°F).

STABILITY

The shelf life of the Rapid fFN Cassette is 18 months from the date of manufacture. Unopened cassettes may be used until the expiration date printed on the foil pouch and the box containing the pouched cassettes. Once the foil pouch is opened, the Rapid fFN Cassette should be used immediately.

MATERIALS PROVIDED

Rapid fFN for the TLiQ® System, REF 01200 (includes Cassettes and Directional Insert).

MATERIALS REQUIRED BUT NOT PROVIDED

1. TLiQ System, REF 01202 (includes Analyzer, Printer, User Manual, and TLiQ QCette®)
2. Rapid fFN Control Kit, REF 01166 (includes Positive Control, Negative Control, and Directional Insert)
3. 200 µL pipette

SPECIMEN COLLECTION

The Hologic Specimen Collection Kit for fetal fibronectin testing is the only acceptable specimen collection system that can be used to collect specimens for this assay. See the Specimen Collection Kit directional insert for complete instructions.

PROCEDURE

PERFORMING ANALYZER QUALITY CONTROL

Use the TLiQ QCette to ensure proper function of the TLiQ Analyzer. *See the TLiQ QCette directional insert for complete instructions.*

SETTING CALIBRATION FOR A RAPID fFN CASSETTE LOT

Select SET CALIBRATION from the TLiQ Analyzer Main Menu and enter the information requested (Cassette Lot # and Calibration Code #). The Cassette Lot # is located on the cassette pouch. The Calibration Code # is located on the cassette box. *See the TLiQ System User Manual for details.*

RUNNING RAPID fFN CONTROL KIT

The Rapid fFN Control Kit must be run each time a new lot or a new shipment of Rapid fFN Cassettes is received. Run the liquid controls as if testing patient samples. *See the Rapid fFN Control Kit directional insert for complete instructions.*

Note: For your convenience, space is provided on the Rapid fFN Cassette Kit box for control testing documentation.

SPECIMEN PREPARATION

Note: Handle the Specimen Transport Tube and all Patient Specimens as if potentially infectious.

1. Allow all Specimen Transport Tubes to come to room temperature before testing.
2. Gently mix the Specimen Transport Tube prior to removing the swab.
3. Open the Specimen Transport Tube cap and swab assembly. The swab shaft should be seated in the cap. Express as much liquid as possible from the swab by rolling the tip against the inside of the tube. Dispose of the used swab in a manner consistent with handling biohazardous materials.

Rapid fFN® Cassette Kit

TESTING PATIENT SAMPLES

Incubation Mode – Internal

Note: The default setting for the TLiQ Analyzer is Internal Incubation Mode. In this mode, the analyzer will time the incubation and automatically initiate reading of the cassette when incubation is complete.

1. Prepare Patient Samples according to the Specimen Preparation section. Mix patient samples before testing.
2. Remove one Rapid fFN Cassette from the foil pouch.
3. Select VIEW SETUP from the TLiQ Analyzer Main Menu to determine if the analyzer is set to Internal Incubation Mode. If Internal Mode is indicated, proceed to step 4. If the analyzer is not set to Internal Mode, select CHANGE SETUP from the Main Menu and change to Internal Incubation Mode. *See the TLiQ System User Manual for details.*
4. Select TEST PATIENT from the TLiQ Analyzer Main Menu and enter the necessary information until the analyzer prompts for cassette insertion.
5. Insert the cassette into the analyzer and press ENTER.
6. When prompted by the analyzer, pipette 200 µL of patient sample into the sample application well of the Rapid fFN Cassette. Immediately press ENTER to activate the analyzer.
7. The analyzer will count down for 20 minutes and analyze the Rapid fFN Cassette.
8. The fFN result for the patient sample will be displayed on the TLiQ Analyzer display screen as POSITIVE, NEGATIVE, or INVALID.
9. If an INVALID result is obtained, retest with 200 µL of additional sample, if available, on a new cassette. If the problem is not corrected, see the TLiQ System User Manual for details, or contact Hologic for technical assistance.

Incubation Mode – External

Note: In External Incubation Mode, the user is responsible for timing the incubation and starting the analysis. If additional cassettes are run, wait 5 minutes before adding sample to the next cassette to allow for the analysis of the previous cassette (approximately 3 minutes) and for entering menu information required by the analyzer for the next cassette.

1. Prepare Patient Samples according to the Specimen Preparation section. Mix patient samples before testing.
2. Remove one Rapid fFN Cassette from the foil pouch.
3. Select VIEW SETUP from the TLiQ Analyzer Main Menu to determine if the analyzer is set to External Incubation Mode. If External Mode is indicated, proceed to step 4. If the analyzer is not set to External Mode, select CHANGE SETUP from the Main Menu and change to External Incubation Mode. *See the TLiQ System User Manual for details.*
4. Select TEST PATIENT from the TLiQ Analyzer Main Menu and enter the necessary information until the analyzer prompts for cassette insertion.
5. Pipette 200 µL of patient sample into the sample application well of the Rapid fFN Cassette and allow to incubate at room temperature for 20 minutes.
6. When the incubation time is complete, insert the cassette into the analyzer and press ENTER. The analyzer will complete the analysis of the Rapid fFN Cassette.
7. The fFN result for the patient sample will be displayed on the TLiQ Analyzer display screen as POSITIVE, NEGATIVE, or INVALID.
8. If an INVALID result is obtained, retest with 200 µL of additional sample, if available, on a new cassette. If the problem is not corrected, see the TLiQ System User Manual for details, or contact Hologic for technical assistance.

INTERPRETATION OF RESULTS

The fFN result for the patient sample will be displayed on the TLiQ Analyzer display screen as POSITIVE, NEGATIVE, or INVALID. The result is positive if the value derived from the patient sample is greater than or equal to the reference calibration value specified by the calibration code. The result is negative if the value derived from the patient sample is less than the reference calibration value specified by the calibration code. The result is invalid if the test does not meet internal quality controls.

Rapid fFN® Cassette Kit

QUALITY CONTROL PROCEDURES

Current Good Laboratory Practice includes the daily use and documentation of either liquid controls or electronic (internal) controls to assure that the calibration of the diagnostic device is maintained within acceptable limits.

The Rapid fFN Control Kit (REF 01166) contains two liquid controls: one Rapid fFN Positive Control and one Rapid fFN Negative Control. The controls are recommended for use in monitoring the performance of the Rapid fFN Cassette. The recommended frequency of use of the controls is once each time a new lot or a new shipment of cassettes is received, or whenever there is uncertainty about Rapid fFN Cassettes. Deviation from the recommended frequency of quality control testing must be validated by the laboratory. If the criteria for controls are not met, do not test patient samples until acceptable results are obtained.

The TLiQ QCette is a quality control device used to verify that the TLiQ Analyzer performs within specification. The TLiQ QCette is a Rapid fFN Cassette replica, containing a membrane with printed test and control lines, which is read by the TLiQ Analyzer. Three different levels of response are measured with this QC device:

1. **High Level:** The blue line at the procedural control position, which is in the high positive range, must be above a minimum threshold value for QC to pass.
2. **Low Level:** The blue printed line at the test line position is in the cutoff range. This line is measured and compared with a value established during instrument setup and must be within 5% of that value for QC to pass.
3. **Negative:** The white space between the blue lines is measured and should always be in the negative range for QC to pass.

Internal controls monitor all components of the TLiQ System and are performed automatically with every test. These internal controls check for: (1) a threshold level of signal at the procedural control position, (2) proper sample flow across the Rapid fFN Cassette, (3) absence of conjugate aggregation (Cassette: Pass/Fail), and (4) proper function of analyzer hardware (Analyzer: Pass/Fail).

LIMITATIONS

The Rapid fFN result should not be interpreted as absolute evidence for the presence or absence of a process that will result in delivery in ≤ 7 days or ≤ 14 days from specimen collection in symptomatic women or delivery in ≤ 34 weeks, 6 days in asymptomatic women evaluated between 22 weeks, 0 days and 30 weeks, 6 days of gestation. A positive fFN result may be observed for patients who have experienced cervical disruption caused by, but not limited to, events such as sexual intercourse, digital cervical examination, or vaginal probe ultrasound. The Rapid fFN result should always be used in conjunction with information available from the clinical evaluation of the patient and other diagnostic procedures such as cervical examination, cervical microbiological culture, assessment of uterine activity, and evaluation of other risk factors.

- **Test results may not be interpreted visually and must be based on the use of the TLiQ System.**
- Modification of the assay protocol described herein may yield erroneous results.
- The assay has been optimized with specimens taken from the posterior fornix of the vagina or the ectocervical region of the external cervical os. Samples obtained from other locations should not be used.
- The safety and effectiveness of using a cutoff other than that provided by the Rapid fFN Cassette Calibration Code has not been established.
- Assay interference from the following components has not been ruled out: douches, white blood cells, red blood cells, bacteria, and bilirubin.
- The presence of infections has not been ruled out as a confounding factor to risk of preterm delivery.
- Assay interference from semen has not been ruled out. Specimens should not be collected less than 24 hours after intercourse. However, even when a patient reports having had intercourse in the previous 24 hours, a negative fetal fibronectin test result is valid.
- Specimens should be obtained prior to digital examination or manipulation of the cervix. Manipulations of the cervix may lead to false positive results.
- Patients with suspected or known placental abruption, placenta previa, or moderate or gross vaginal bleeding should not be tested.

Rapid fFN® Cassette Kit

EXPECTED VALUES

Among symptomatic women, elevated levels ($\geq 0.050 \mu\text{g/mL}$) of fFN between 24 weeks, 0 days and 34 weeks, 6 days indicate increased risk of delivery in ≤ 7 or ≤ 14 days from sample collection. Similarly, among asymptomatic women, elevated levels of fFN between 22 weeks, 0 days and 30 weeks, 6 days indicate increased risk of delivery in ≤ 34 weeks, 6 days of gestation. The cutoff of $0.050 \mu\text{g/mL}$ fFN was established in a multicenter study conducted to evaluate the association between fetal fibronectin expression during pregnancy and preterm delivery (5).

Only subjects with symptoms of preterm labor or preterm rupture of membranes were eligible for this study. Of the total study population, the association of vaginal fetal fibronectin expression to preterm delivery was evaluated for 117 symptomatic women with intact amniotic membranes. The strength of this association was determined at a variety of cutoffs using receiver operator characteristic (ROC) curves. These results show that the optimal sensitivity and specificity is provided at a cutoff of $0.050 \mu\text{g/mL}$ fFN. Subsequent studies evaluating the fFN Enzyme Immunoassay test as a predictor of preterm delivery among symptomatic women and asymptomatic women with singleton pregnancies have confirmed that the optimal cutoff is $0.050 \mu\text{g/mL}$ fFN (5,13,15,18). Laboratory studies have confirmed the optimal cutoff of $0.050 \mu\text{g/mL}$ fFN in the Rapid fFN. This analytical cutoff was used to show equivalency between the fFN Enzyme Immunoassay and the Rapid fFN.

PERFORMANCE CHARACTERISTICS

Accuracy

A comparison of the Fetal Fibronectin Enzyme Immunoassay to the Rapid fFN was assessed in 587 cervicovaginal samples. The two test formats were equivalent 94.9% of the time. The Kappa coefficient for inter-test agreement was 0.81 with 95% confidence interval of [0.75, 0.88].

Table
Analytical Concordance Between fFN Enzyme Immunoassay and Rapid fFN (n=587)
(Symptomatic and Asymptomatic Women)

	fFN Enzyme Immunoassay (+)	fFN Enzyme Immunoassay (-)	
Rapid fFN (+)	77	17	94
Rapid fFN (-)	13	480	493
	90	497	587

Within-Run Reproducibility

Within-run reproducibility was determined using three lots of Rapid fFN Cassettes. Twenty replicates each of the Rapid fFN Negative Control (fFN concentration approximately $0.015 \mu\text{g/mL}$) and of the Rapid fFN Positive Control (fFN concentration approximately $0.080 \mu\text{g/mL}$) were tested on three lots of Rapid fFN Cassettes using three different TLiQ Analyzers. (One lot of cassettes was run on one analyzer on one day, on three different occasions.) The results from these studies show that the Rapid fFN correctly identified all specimens 100% of the time.

Between-Run Reproducibility

Between-run reproducibility was determined in 36 independent assays, using different lot combinations of Rapid fFN Cassettes (three lots) and TLiQ Analyzers (six each). Each assay consisted of the Rapid fFN Negative Control and assay of the Rapid fFN Positive Control. The Rapid fFN correctly identified all samples 100% of the time.

Interfering Substances

Care must be taken not to contaminate the swab or cervicovaginal secretions with lubricants, soaps, disinfectants, or creams (e.g., K-Y® Jelly lubricant, Betadine® disinfectant, Monistat® cream). Lubricants or creams may physically interfere with absorption of the specimen onto the swab. Soaps or disinfectants may interfere with the antibody-antigen reaction.

Various concentrations of pharmacologic agents were added to specimens containing approximately $0.015 \mu\text{g/mL}$ to $0.080 \mu\text{g/mL}$ fFN and assayed in triplicate. The drugs added were: **prostaglandin E₂** (up to $250 \mu\text{g/mL}$), **ampicillin** (up to $100 \mu\text{g/mL}$), **cephalexin** (up to $18 \mu\text{g/mL}$), **erythromycin** (up to $10 \mu\text{g/mL}$), **gentamycin** (up to $4 \mu\text{g/mL}$), **dexamethasone** (up to $200 \mu\text{g/mL}$), **magnesium sulfate** (up to $50 \mu\text{g/mL}$), **oxytocin** (up to $100 \mu\text{g/mL}$), **terbutaline** (up to $100 \mu\text{g/mL}$), and **ritodrine** (up to $10 \mu\text{g/mL}$). These drugs did not interfere with the assay at the concentration limits cited above.



Rapid fFN® Cassette Kit

BIBLIOGRAPHY

1. American College of Obstetricians and Gynecologists. Preterm Labor. **Technical Bulletin, Number 133**, October, 1989.
2. Creasy RK, Resnick R. **Maternal and Fetal Medicine: Principles and Practice**. Philadelphia: W.B. Saunders; 1989.
3. Creasy RK, Merkatz IR. Prevention of preterm birth: clinical opinion. **Obstet Gynecol** 1990;76(Suppl 1):2S–4S.
4. Morrison JC. Preterm birth: a puzzle worth solving. **Obstet Gynecol** 1990;76(Suppl 1):5S–12S.
5. Lockwood CJ, Senyei AE, Dische MR, Casal DC, et al. Fetal fibronectin in cervical and vaginal secretions as a predictor of preterm delivery. **New Engl J Med** 1991;325:669–74.
6. Maymon R, Bahari C, Moroz C. Placental isoformin measured by a specific monoclonal antibody as a predictive marker for preterm contraction outcome. **Obstet Gynecol** 1989;74:597–9.
7. Wasmoen TL, Coulam CB, Leiferman KM, Gleich GJ. Increases of plasma eosinophil major basic protein levels late in pregnancy predict onset of labor. **Proc Natl Acad Sci USA** 1987;84:3029–32.
8. Matsuura H, Hakomori SI. The oncofetal domain of fibronectin defined by the monoclonal antibody FDC-6: its presence in fibronectins from fetal and tumor tissues and its absence in those from normal adult tissues and plasma. **Proc Natl Acad Sci USA** 1985;82:6517–21.
9. Matsuura H, Takio K, Titani K, Greene T, et al. The oncofetal structure of human fibronectin defined by monoclonal antibody FDC-6. Unique structural requirement for the antigen specificity provided by a glycosylhexapeptide. **J Biol Chem** 1988;263:3314–22.
10. Feinberg RF, Kliman HJ, Lockwood CJ. Is oncofetal fibronectin a trophoblast glue for human implantation? **Am J Pathol** 1991;138:537–43.
11. Morrison JC, Allbert JR, McLaughlin BN, Whitworth NS, et al. Oncofetal fibronectin in patients with false labor as a predictor of preterm delivery. **Am J Obstet Gynecol** 1993;168:538–42.
12. Inglis SR, Jeremias J, Kuno K, Lescale K, et al. Detection of tumor necrosis factor- α , interleukin-6, and fetal fibronectin in the lower genital tract during pregnancy: relation to outcome. **Am J Obstet Gynecol** 1994;171:5–10.
13. Iams J, Casal DC, Goodwin TM, Kreaden US, et al. Fetal fibronectin improves the accuracy of diagnosis of preterm delivery. **Am J Obstet Gynecol** 1995;173:141–5.
14. Burrus DR, Ernest JM, Veille JC. Fetal fibronectin, interleukin-6, and C-reactive protein are useful in establishing prognostic sub-categories of idiopathic preterm labor. **Am J Obstet Gynecol** 1995;173:1258–62.
15. Bartnicki J, Casal DC, Kreaden US, Saling E, Vetter K. Fetal fibronectin in vaginal specimens predicts preterm delivery and very low birth weight infants. **Am J Obstet Gynecol** 1996;174:971–4.
16. Lockwood CJ, Wein R, Lapinski R, Casal D, et al. The presence of cervical and vaginal fetal fibronectin predicts preterm delivery in an inner-city obstetric population. **Am J Obstet Gynecol** 1993;169:798–804.
17. Leeson SC, Maresh MJA, Martindale EA, Mahmood T, et al. Detection of fetal fibronectin as a predictor of preterm delivery in high risk asymptomatic pregnancies. **Br J Obstet Gynecol** 1996;103:48–53.
18. Goldenberg RL, Mercer BM, Meis PJ, Copper RL, et al. The preterm prediction study: fetal fibronectin testing and spontaneous preterm birth. **Obstet Gynecol** 1996;87:643–8.
19. Morrison JC, Naef RW, Botti JJ, Katz M, et al. Prediction of spontaneous preterm birth by fetal fibronectin and uterine activity. **Obstet Gynecol** 1996;87:649–55.
20. McGregor JA. Salivary estriol as risk assessment for preterm labor: a prospective trial. **Am J Obstet Gynecol** 1995;173:1337–42.

This product may be covered by one or more U.S. patents identified at <http://hologic.com/patentinformation>

© 2017 Hologic, Inc. All rights reserved.

Hologic, Rapid fFN, QCette and/or associated logos are trademarks and/or registered trademarks of Hologic, Inc. and/or its subsidiaries in the United States and/or other countries. All other trademarks, registered trademarks, and product names are the property of their respective owners.

Rapid fFN® Cassette Kit

TECHNICAL SERVICE AND ORDERING INFORMATION

USA/CANADA ONLY

Tel: 1-800-442-9892






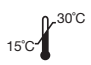



Fax: 1-508-263-2956

ALL OTHER COUNTRIES

Tel: +1-508-263-2900

For additional contact information, go to www.ffntest.com

Following is an explanation of the symbols that may appear on your product.

Symbol	Title	Description	Standard information
	Caution, consult instructions for use	Indicates the need for the user to consult the instructions for use for important cautionary information such as warnings and precautions that cannot, for a variety of reasons, be presented on the medical device itself.	ISO 15223-1 Medical devices—Symbols to be used with medical device labeling and information to be supplied, Section 5.4.4
	Consult instructions for use	Indicates the need for the user to consult the instructions for use.	ISO 15223-1 Medical devices—Symbols to be used with medical device labeling and information to be supplied, Section 5.4.3
	Use-by date	Indicates the date after which the medical device is not to be used.	ISO 15223-1 Medical devices—Symbols to be used with medical device labeling and information to be supplied, Section 5.1.4
	<i>In vitro</i> diagnostic medical device	Indicates a medical device that is intended to be used as an <i>in vitro</i> diagnostic medical device	ISO 15223-1 Medical devices—Symbols to be used with medical device labeling and information to be supplied, Section 5.5.1
	Batch code	Indicates the manufacturer's batch code so that the batch or lot can be identified	ISO 15223-1 Medical devices—Symbols to be used with medical device labeling and information to be supplied, Section 5.1.5
	Temperature limit	Indicates the upper and lower limit of temperature to which the medical device can be safely exposed.	ISO 15223-1 Medical devices—Symbols to be used with medical device labeling and information to be supplied, Section 5.3.7
	Catalogue Number	Indicates the manufacturer's catalogue number so that the medical device can be identified	ISO 15223-1 Medical devices—Symbols to be used with medical device labeling and information to be supplied, Section 5.1.6
	Do not re-use	Indicates a medical device that is intended for one use, or for use on a single patient during a single procedure.	ISO 15223-1 Medical devices—Symbols to be used with medical device labeling and information to be supplied, Section 5.4.2
	Manufacturer	Indicates the medical device manufacturer, as defined in the EU Directives 90/385/EEC, 93/42/EEC and 98/79/EC.	ISO 15223-1 Medical devices—Symbols to be used with medical device labeling and information to be supplied, Section 5.1.1